


BABY BLUE

BRINGING THE SPA HOME



USER MANUAL

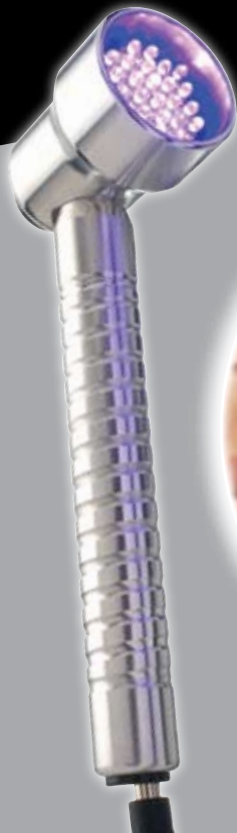


MAKING ACNE SING THE BLUES

REFINE YOUR SKIN...ON YOUR TIME.

DESIGNED SPECIFICALLY TO CLEAR AND CALM IRRITATED AND STRESSED SKIN, BABY BLUE IS THE MOST POWERFUL ANTI-ACNE SKINCARE DEVICE AVAILABLE FOR AT-HOME USE.

ACNE IS CAUSED BY THE BODY PRODUCING TOO MUCH SEBUM, THE SLIGHTLY OILY SUBSTANCE PRODUCED BY THE SEBACEOUS GLANDS IN HAIR FOLLICLES. THE BACTERIA GET INTO THE FOLLICLES WHERE IT FLOURISHES AND CREATES INFECTION. THIS CAN TURN INTO A WHITEHEAD, BLACKHEAD OR BLEMISH. INFLAMMATION OFTEN FOLLOWS CREATING A RED, UNSIGHTLY APPEARANCE.

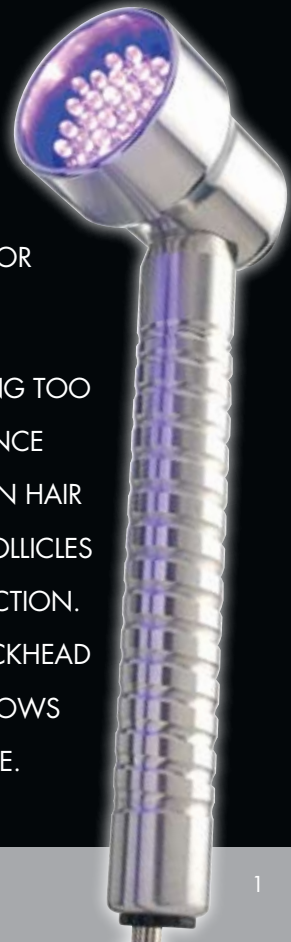

BABY BLUE

BABY BLUE



THE BABY BLUE LIGHT WORKS by killing the acne-causing bacteria, *P. acnes*, and treating inflammatory acne vulgaris. The *P. acnes* bacteria releases small molecules called porphyrins. When porphyrins are exposed to the baby blue's wavelengths of natural light, they produce free radicals that kill the bacteria.

THE BABY BLUE LIGHT:

- Reduces blemishes
- Kills acne-causing bacteria (*Propionibacterium acnes*)

While acne is not considered a serious health threat, severe acne can result in permanent scarring and can be a source of significant emotional distress.

DELIVERING RESULTS FROM THE INSIDE OUT, Baby Blue is 100% natural and requires no harsh chemicals or medication. It is safe, effective, painless, non-invasive and requires no downtime. It is so gentle, it can be used on all skin types and skin tones, including sensitive and irritated skin.

PREPARATION

A clean face is important in Light Therapy treatment. Remove all makeup and jewelry from areas to be treated. Clean the area to be treated with your favorite pH balanced cleanser.

WARNING

- Do not scratch or damage the light (LEDs)
- Keep away from flammable aerosol products being used or where oxygen is being administered.
- Keep the cord away from water and heat.
- Do not operate this device if the cord or the plug has been damaged or submerged in water
- Do not put into dishwasher, submerge in water or cleaning solutions

CLEANING

The LEDs should be cleaned following each use. After unplugging the unit, use a soft cloth, and clean with rubbing alcohol. You may use cotton swabs to gently clean around the LEDs.

PACKAGE CONTENTS



Baby blue

- Travel pouch
- User manual
- International power supply

Medically proven wavelength of 415 made of solid aluminum. |

TREATMENT RECOMMENDATION

Use the Baby Blue up to 15 minutes per day on acne prone areas. You may use the baby blue longer for problem areas.

(Unlike the Baby Quasar, the Baby Blue does not have a soft tone timer.) Use the Baby Blue on days you're not treating with the Baby Quasar (red) for up to 15 minutes per day on acne prone areas. Combining the Baby Blue with the Baby Quasar is the most powerful acne treatment. Use the Baby Quasar (red) on the continuous wave mode. (green status indicator light) for four minutes in each section as outlined in the Baby Quasar



BABY BLUE

BABY QUASAR

(pages 12–13).

Light treatment with the Baby Blue is recommended to be used 2 to 3 days per week for approximately 15 minutes per day.

THE BABY BLUE should be used once or twice weekly for 15 minutes to control acne bacteria. With regular and continued use, baby blue's light therapy treatment will dramatically improve the appearance of acneic skin, resulting in a clean, clear complexion. For best results, combine with the baby quasar's red and infrared light to help control sebum, stimulate healing, and reduce inflammation. Dermatologist recommended.

WARRANTY

Quasar Light Therapy warrants this product to be free from defective workmanship and materials, and agrees that it will, at its option, either repair or replace the defective product or part thereof with a new or remanufactured equivalent at no charge to

QUASAR LIGHT THERAPY

1431 Tallevast Rd.

Sarasota, FL 34243

(941) 306-5812

www.babyquasar.com

BABY QUASAR

BRINGING THE SPA HOME

COME TO THE LIGHT.

GORGEOUS, HEALTHIER, YOUNGER-LOOKING
SKIN IS ONLY A FEW MINUTES A DAY AWAY.

Using superior technology, baby quasar works by stimulating the production of collagen, thereby diminishing the appearance of fine lines, wrinkles and enlarged pores.

